


Please read the Instructions carefully and consult your care provider if necessary before continuing.
READ THE BACK FOR CAUTIONS AND WARNINGS

SETTING UP THE MEDCENTER SYSTEM:

- Remove each individual Pill Cassette (dated “1” to “31”) and lay them out with lids open in numerical order.
- Follow each pill’s prescribed instructions and place each pill into the proper compartment for each day.
- The Pill Cassettes should be inserted into the MedCenter with the green end up denoting they are loaded and ready.
- Pull the plastic tab on the back of the MedCenter Clock to connect the batteries and activate the clock. At activation, **the clock has 4 default alarm times (9:00am, 1:00pm, 6:00pm, 10:00pm)**, if necessary, set the alarms to your desired times (see SETTING THE ALARMS) and place Clock in the MedCenter Stand as shown.

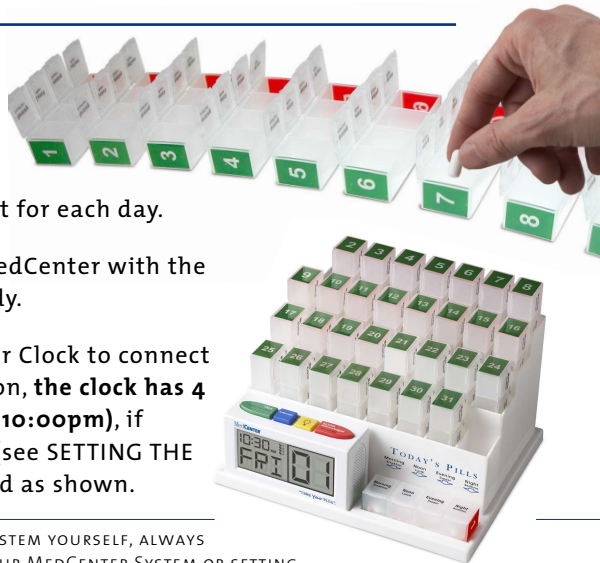

NOTE: IN CASE YOU ARE NOT CAPABLE OF FILLING YOUR MEDCENTER SYSTEM YOURSELF, ALWAYS HAVE YOUR CARE PROVIDER OR PHARMACIST ASSIST YOU IN FILLING YOUR MEDCENTER SYSTEM OR SETTING YOUR MEDCENTER CLOCK ACCORDING TO YOUR DOCTOR OR HEALTH CARE PROVIDER’S INSTRUCTIONS EACH MONTH.

 **CAUTION: *Please Read Carefully***

The MedCenter System is **to be used only as an organizational aid in taking your medications**. The responsibility for correct medication compliance is solely up to the consumer. The makers of MedCenter Systems are not responsible for incorrect loading or unintended use of the pill cassette system, for improper setting or battery life on the MedCenter Clock or its alarm functions, for incorrect consumption of medication due to any use of the MedCenter System.

Carefully **check prescribed instructions for each medication BEFORE organizing your pills in the MedCenter System**. If you are not capable of filling your MedCenter System yourself as prescribed, always have a care provider or pharmacist assist you.

Store the MedCenter Clock and MedCenter System away from heat and water.

In order to maintain proper operation, Batteries should be replaced once every year. The consumer is responsible for the replacement of the battery and setting of the clock.

MedCenter Talking Clock will ONLY remind you if you can hear the alarm. Always maintain a backup method to remind you to take your pills.

 **WARNING: *Please Read Carefully***

Keep the MedCenter Clock and MedCenter System OUT OF REACH OF CHILDREN.

Always verify the medication type, dosage and schedule with Doctor, Pharmacist or Care Provider BEFORE taking any medication.